



# CombiCarrier II<sup>®</sup>

## Scoop Stretcher & Extrication Board

**REVOLUTIONARY**  
Patient Handling



The CombiCarrierII<sup>®</sup> functions as a scoop-type stretcher minimizing patient movement with no logrolling required and also serves as an improved extrication board during vehicle extrication situations. The CombiCarrierII meets the latest Spinal Motion Restriction protocols providing optimal patient handling in the most extreme environments including water rescue.

The concave patient surface limits lateral movement of the victim providing more secure stabilization. The split-litter design also allows for easy removal once the patient is on the ambulance cot or delivered to the emergency department. Built to take a beating, many of the original units are still in operation with zero failures in over 15 years.

The innovative design of the CombiCarrierII provides emergency medical and rescue personnel with two products in one, saving valuable storage space, weight and critical budget dollars. **Comes with an industry leading 5-year warranty.**



*Safely extricate and transfer a patient to the cot. The CombiCarrierII is easily removed eliminating the loss of your valuable equipment.*



### QUICK AND EASY

Simple intuitive design minimizes training of all EMS, hospital and athletic training personnel.



### AUTO-LOCK LATCH SYSTEM

Off-set sturdy latch is built for reliability and durability making it much easier to apply and remove the device under patient load.



### VERSATILE PATIENT TRANSFER

Can be separated at either end or split into two halves, eliminating unnecessary patient movement and making it ideal for patients with suspected hip, pelvic and spinal injuries.



### CONCAVE PATIENT SURFACE

Provides excellent lateral support with neutral alignment. Minimizes movement of patients spine during lifting and transferring. Continuous inner edge eliminates pinching.



### SEAMLESS PLASTIC CONSTRUCTION

No rivets, drive screws or roll pins. High density polyethylene is easily cleaned and disinfected. Helps comply with CDC and OSHA guidelines.




### X-RAY COMPATIBLE

X-ray compatible. Easy to remove for highly detailed diagnostic imaging requirements.

**HARTWELL MEDICAL**

# CombiCarrier II<sup>®</sup>

## Scoop Stretcher & Extrication Board

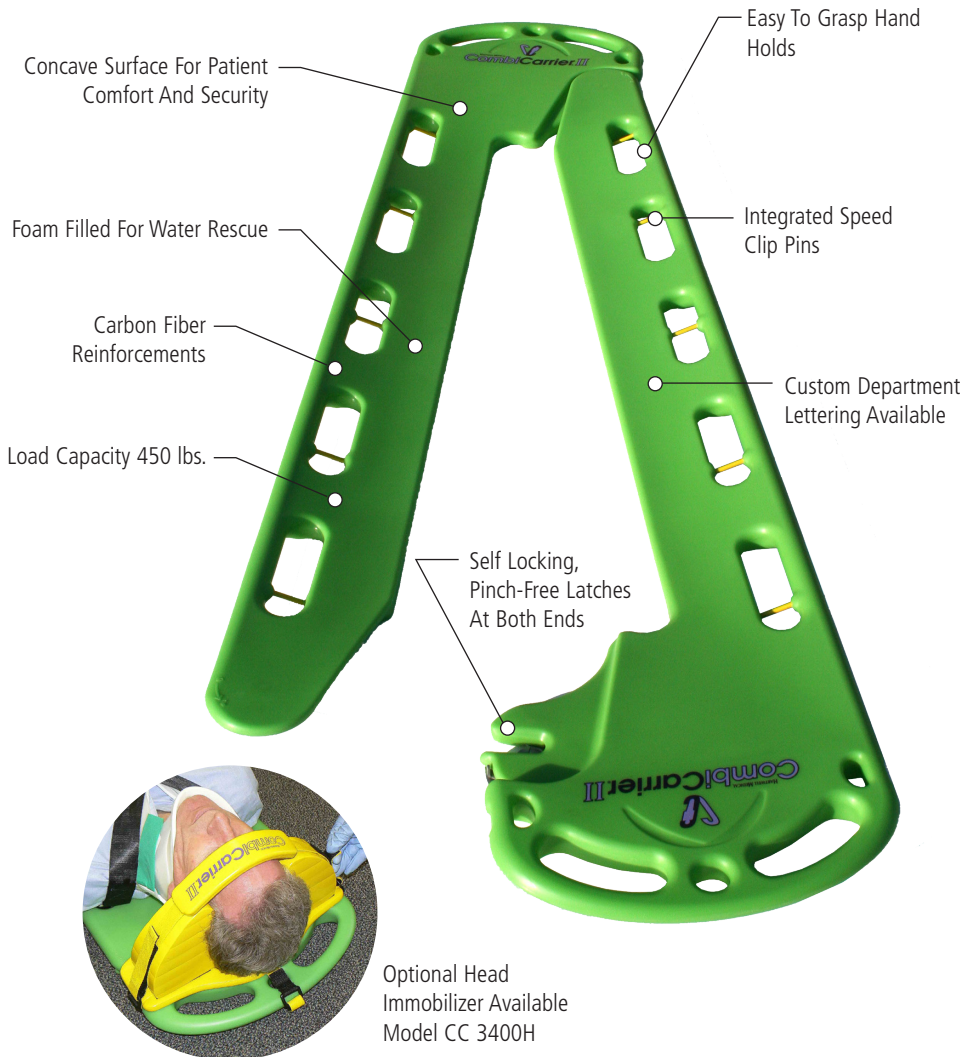

Functions Like A Scoop Stretcher



Tapered Ends To Simplify Extrication



Fits In A Standard Backboard Compartment



Can Be Easily Separated At Either End

### SPECIFICATIONS - Model CC 2200P and CC 2200P-OD

Weight	15.6 lbs. (7.1 kg)
Dimensions	73.50"L x 16.50"W x 2.19"D (186.7 cm x 41.9 cm x 5.6 cm D)
Load Capacity	450 lbs. (205 kg)
Shell	High density polyethylene
Pins	Fiberglass
Filling	Urethane foam
Latches	Aluminum, stainless steel

For optional ordering configurations please see our website. Tactical version available in olive drab.



Contact your local dealer or visit us online at [www.HartwellMedical.com](http://www.HartwellMedical.com)  
For detailed product information, training and education materials please visit our website