


S-Guide

Malleable Intubating Guide

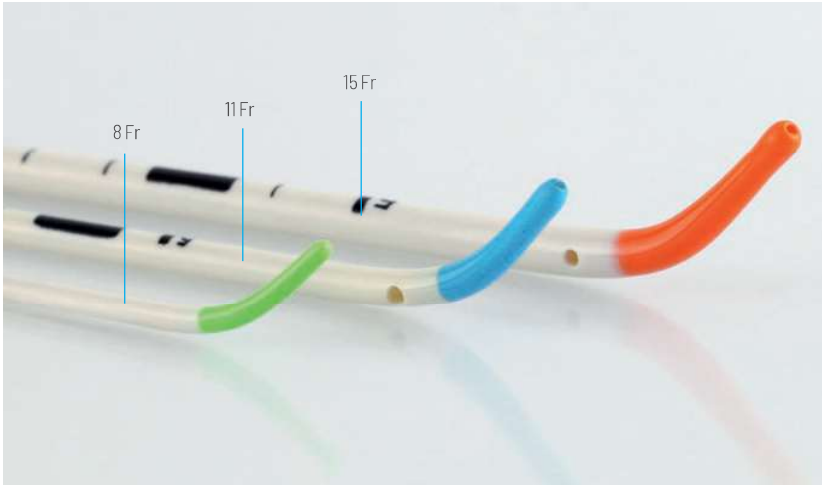


The S-Guide is an intubation aid for difficult airway management and facilitates orotracheal intubation, especially in patients where visualisation of the glottis is inadequate.

Due to its **malleability**, it is suitable for direct and indirect laryngoscopy, in particular for video laryngoscopes with non channelled blades.

Supplied with two types of O₂ Connectors, **oxygenation*** is possible throughout the whole intubating process, not impairing oxygenation even when the S-Guide is bent. For patient safety the S-Guide is provided with a coloured, pre-formed and atraumatic **soft tip**.

*Oxygenation is possible with the sizes 11 Fr and 15 Fr. Oxygenation is not possible with the size 8 Fr.


Features

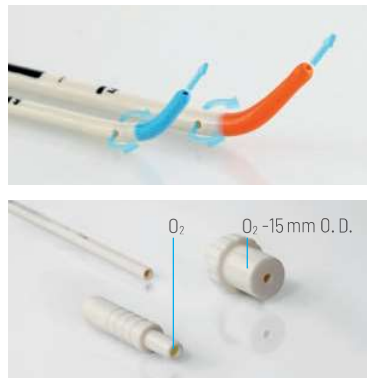

Malleability

Oxygenation*

Soft tip

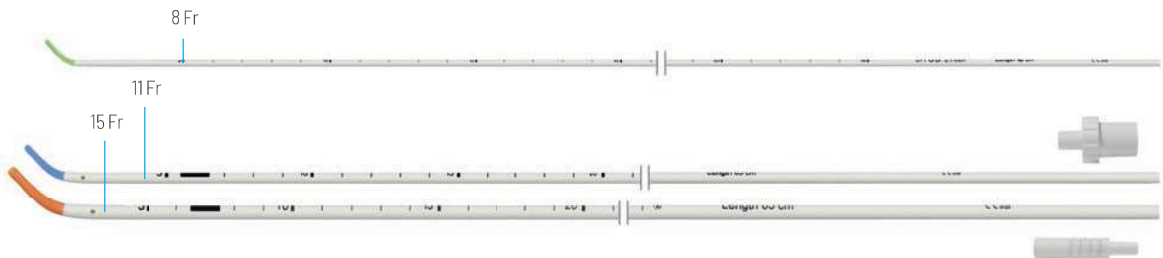

Order information

S-Guide / For single use, sterile
 Centimetre marking from 5 to 40 cm with the sizes 11 Fr and 15 Fr,
 two types of O₂ Connectors included
 Centimetre marking from 5 to 30 cm with the size 8 Fr

Size	Length	For Tracheal Tube	REF	Box
8 Fr	40 cm	≥ 3.0 mm I.D.	33-92-400-1	5
11 Fr	65 cm	≥ 4.5 mm I.D.	33-91-650-1	5
15 Fr	65 cm	≥ 6.0 mm I.D.	33-90-650-1	5

Additional information



Application Video



Brochure



Recommendation of use

References see page 28