

Contour and straight shape

Colour coding

Low profile design



Tourniquet Dispo Cuff

The hygienic single use solution

Tourniquet Dispo Cuff


The Tourniquet Dispo Cuff is a sterile, single use product especially suitable for contaminated procedures, preventing against HAI (Hospital-acquired infection). The Dispo Cuff is available in straight shape (for small or midsize upper extremity) and contour shape (mainly for lower conical shaped extremity) and is provided in different sizes, which are colour coded for optimal cuff selection.

Features



Cuff-Design

- Large hook and loop fastener for reliable cuff fixation
- Tubing made of kink-resistant silicone material
- Low profile design and light weight
- Contour and straight shape
- PVC-free

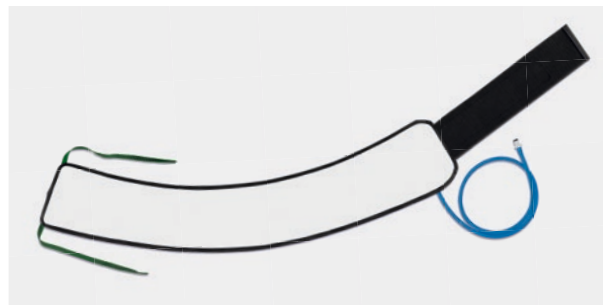


Contour shape

- Conical design with a smaller diameter distally than proximally when applied
- Occludes the blood flow at lower cuff pressures due to an improved fit to conical extremities
- Integrated stiffener to prevent cuff migration



Packaging

- Packed in dispenser box
- Colour coding makes size selection and ordering easier
- Easy storage due to low profile packaging

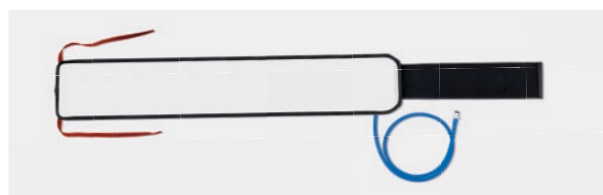

Straight shape

- Cylindrical design
- Does not show any difference in diameter between proximal and distal part

Order information

Dispo Cuff / Single Cuff, for single use, sterile

Length	Colour	Single Port	Dual Port	Box
20 cm	-	REF 20-34-700SLZ-1	REF 20-36-700SLZ-1	10
30 cm	■	REF 20-34-710SLZ-1	REF 20-36-710SLZ-1	10
35 cm	■	REF 20-34-711SLZ-1	REF 20-36-711SLZ-1	10
46 cm	■	REF 20-34-712SLZ-1	REF 20-36-712SLZ-1	10
46 cm	■	REF 20-34-715SLZ-1*	REF 20-36-715SLZ-1*	10
61 cm	■	REF 20-34-722SLZ-1*	REF 20-36-722SLZ-1*	10
76 cm	■	REF 20-34-727SLZ-1*	REF 20-36-727SLZ-1*	10
86 cm	■	REF 20-34-728SLZ-1*	REF 20-36-728SLZ-1*	10
107 cm	■	REF 20-34-729SLZ-1*	REF 20-36-729SLZ-1*	10

Dispo Cuff / Double Cuff for Intravenous regional anaesthesia (IVRA), for single use, sterile

Length	Colour	REF	Box
30 cm	■	20-30-710SLZ-1	10
46 cm	■	20-30-712SLZ-1	10
61 cm	■	20-30-722SLZ-1	10

Extension Tubing / For single use, sterile

Length	REF	Box
100 cm	20-39-100-1	5

*Contour shape

The medical devices in this advertising material are manufactured without the use of natural rubber latex, unless otherwise specified. The medical devices in this advertising material do not contain phthalates which require labelling according to CLP Regulation (EC) 1272/2008.